


# SOFT POINT® DELUXE LEATHER ORTHOTIC INSOLES

**Model No: F35**

**Sizes: B - F**

**Color: Brown**

The Deluxe leather Orthotic Insoles combine the comfort of premium leather and high-density cushioning foam with semi-rigid arch support and motion control. Unique tri-layer design has many features to help relieve and prevent common foot problems. Discover the many benefits for yourself.

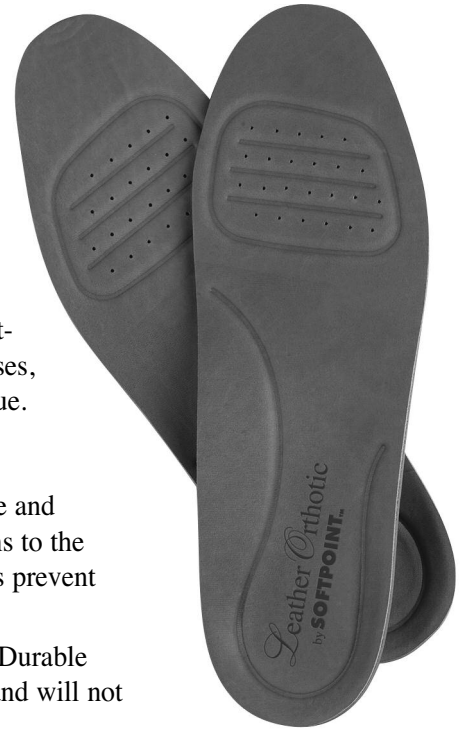

**Indicated for:** Plantar or Arch Pain, Shin Splints, Heel Pain/Spurs, Arch Support, Post-Operative Conditions, Forefoot Pain and Calluses, and Lower Back, Leg, and Foot Pain and Fatigue.

## FEATURES:

- **Premium Leather Top Cover** - Comfortable and breathable. Premium, supple leather conforms to the foot with use for even greater comfort. Helps prevent and relieve calluses on the forefoot and toes.
- **High-Density Closed Cell Cushion Foam** - Durable foam backing is resilient to long term wear and will not break down.
- **Semi-Rigid Plastic Arch Support** - Provides arch support to help alleviate heel and arch pain and forefoot stress. Semi-rigid support also assists in controlling excessive motion and lower leg rotation to reduce the stress on knee and hip joints.
- **Movable Metatarsal Pad** - Includes an optional movable metatarsal pad which can be positioned exactly where needed for additional cushioning under the metatarsals. Helps relieve pain from neuromas or metatarsalgia.
- **Gel Heel Pad** - A gel insert is located at the heel to help reduce stress and cushion heel strikes at the most sensitive area of the foot.

## HOW TO APPLY:

Remove the existing insole from the shoe. Place the Soft Point Insole in the respective shoe with the leather side up. The insole should lay flat if it curls at the front or the sides slight trimming with scissors may be necessary. Only slight trimming from the toe and sides is recommended – do not trim the heel. Use the removed insole and the trim lines as guides to ensure sized correctly. The insoles can be easily moved from one shoe to another. .

## CARE INSTRUCTIONS:

Wipe clean. Air dry.

**SIZING:** Size according to US shoe size.

Size	Men's Shoe	Women's Shoe
B	—	5 - 7
C	5 - 7	7 - 9
D	7 - 9	9 - 11
E	9 - 11	—
F	11 - 13	—

**800-327-4110**



Website: [www.flaorthopedics.com](http://www.flaorthopedics.com)

©2004 FLA ORTHOPEDICS, INC.